



Captair Pyramid

Portable glove bag

Safe and convenient transport, handling, and testing of hazardous samples



Ready-to-use:

- Transport and open suspicious packages
- Sample and user protection
- Protection from potential leaks, or open vials
- Versatility to be used as test chamber

Portable glove bag

Designed to offer protection from hazardous, known, or potentially contaminated samples. Specimens that are obtained can be placed in the Pyramid, sealed and transported to the lab for testing. If in-field testing is required the Pyramid can be used as a glove box and injected with an inert gas for sample stability.

Easy to set up and use in the field or lab

It is a multi-function disposable glove bag made of high-quality transparent PVC. Assembled in a few seconds, it can be used anywhere. The slanted shape of the enclosure provides an ergonomic working position for the user.

Light, mobile and disposable, it is an ideal protection tool which can suit many requirements in the laboratory or in the field.


Optimal air-tightness

Each Pyramid is factory-tested to ensure an airtight seal by inflating with compressed air at a pressure of 2.5 mm Hg.

Once tested, it is given a serial number and a certificate to ensure safety.


Easy transport

Folds into an easy-to-carry size

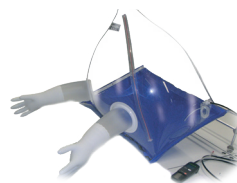

Each unit individually tested and issued a unique serial number for full traceability

Specifications

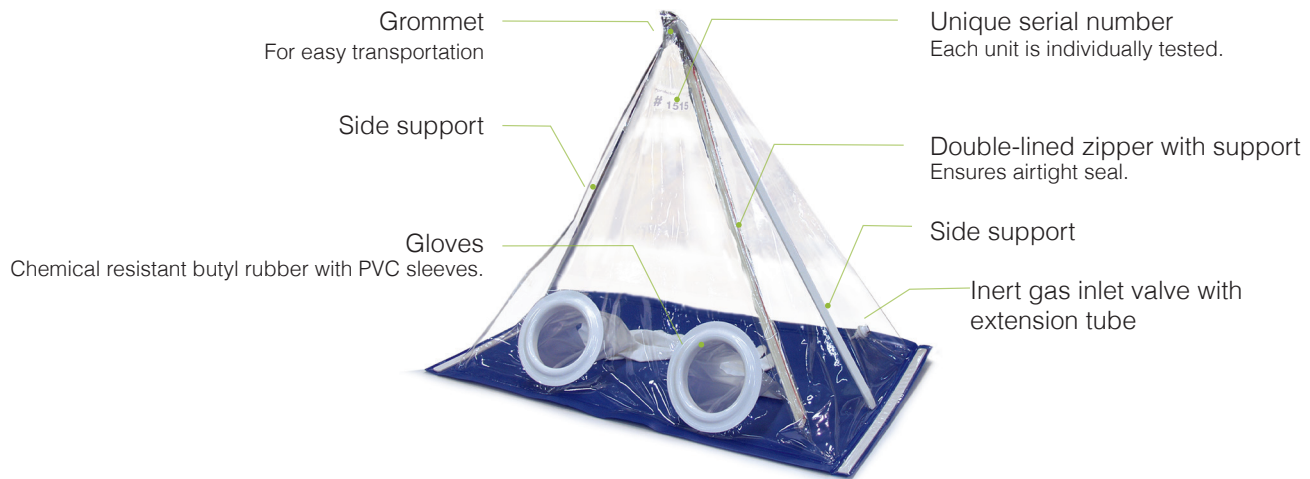

Dimensions

| Width (mm/inches) | Depth (mm/inches) | Height (mm/inches) |
|-------------------------|-------------------|-------------------------|
| 860 / 33 ^{7/8} | 560 / 22 | 725 / 28 ^{1/2} |

Technical specifications

| | |
|------------------------------|---|
| Enclosure and base | Flexible PVC |
| Fastener | Double lip zipper (600 mm / 23 ^{5/8} in) |
| 2 arm openings | Rigid PVC (115 mm / 4 ^{1/2} in Ø) |
| Grommet | 16 mm / 5/8 in Ø |
| Operating temperature | -25°C +45°C / -13°F +113°F |
| Weight | 1.55 kg / 3.42 lbs |

Features

| | |
|----------------------|--|
| Gloves | Chemical resistant butyl rubber with PVC sleeves |
| Serial number | Unique reference number for better traceability |
| Valve | Gas insertion valve |

France
+33 (0) 2 32 09 55 80 | ventes@erlab.net

Germany
0800 330 47 31 | verkauf@erlab.net

United States
+1 800-964-4434 | captairsales@erlab.com

United Kingdom
+44 (0) 1722 341 940 | salesuk@erlab.net

China
+86 (0) 512 5781 4085 | sales.china@erlab.com.cn

Italy
+39 (0) 2 89 00 771 | vendite@erlab.net

Malaysia
+60 (0) 7 3 555 724 | erlab@tm.net.my

Spain
+34 93 673 24 74 | ventas@erlab.net

